**Description:**

The Model 4902 CRYOTherm Single-Tank Console provides comfortable and convenient of cold hydrotherapy. Water in the tank is precisely maintained at selected temperatures with a hermetic refrigeration system and electronic temperature control.

**Specifications:**

Weight / Tanks Empty	2,250 Pounds
Weight / Tanks Full	8,078 Pounds
Tank Dimensions	60" L x 84" W x 34" D
Tank Capacity	698 Gallons
Temperature Range	40° F to 70° F
Tank Construction	12 Gauge, Type 304 Stainless Steel
Exterior Construction	Dupont Cameo White Corian
Voltage	208-240/120 Volt, Single Phase, 60 Hertz
Full Load Amperage	23 Amps
Circuit Breaker	2 Pole, 30 Amp, GFI
Refrigeration Capacity	28,500 BTU @ 50°F
Pump	1.5 HP, Dual Speed
Hydrotherapy Jets	Six per Tank

***Please read thoroughly before use and retain for later reference***



Grimm Scientific Industries, Inc.  
P.O. Box 2143, Marietta, OH 45750  
www.grimmscientific.com

Tel: 800-223-5395  
Fax: 740-374-5745  
e-mail: grimm@ee.net

**Important Safeguards:**


*The following safety precautions should always be followed when using CRYOTherm.*

1. Use CRYOTherm only as directed in this manual. Improper use can cause injury.
2. Do not plug CRYOTherm into any electrical outlet that is not protected by a ground fault interrupter and matches the voltage and amperage specification.
3. Do not modify the power cord in any way or use with any type extension cord.
4. Maintain CRYOTherm as described in this manual. Improper maintenance can cause injury.
5. Use only Grimm-approved parts or service. Improper parts or service can cause injury.
6. Do not utilize warm or hot water treatment where there is hemorrhaging or swelling.
7. Advise patients to refrain from adjusting controls, and never leave patients alone.
8. Thoroughly clean the tanks with a disinfectant between treatments, particularly after treating open lesions (see "Maintenance" section of this manual).
9. An antibacterial agent added when filling the tanks helps prevent or reduce infection.
10. Keep loose material (bandages, gauze, washcloths, etc.) out of tanks as they can block pump intakes.


**Water Temperature Terminology:**

Very Cold	40° F to 55°F	5° C to 13° C
Cold	55° F to 60°F	13° C to 16° C
Warm	90° F to 98°F	32° C to 37° C
Hot	98° F to 104°F	37° C to 40° C
Very Hot	104° F to 110°F	40° C to 43° C

## Control Functions:

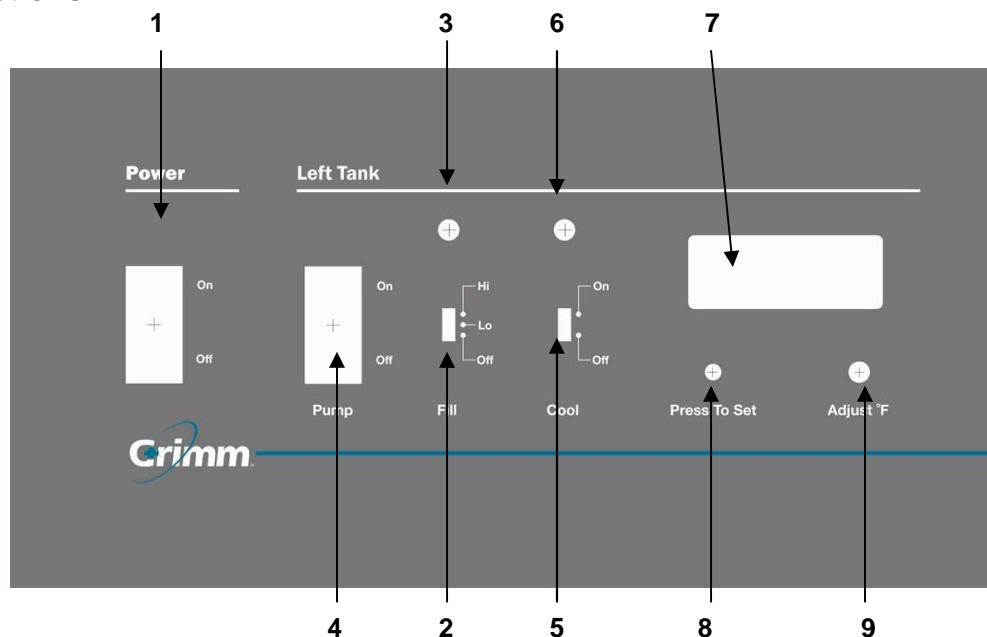

1. **Power Switch** - Activates all electrical components.
2. **Fill Switch** - Activates automatic water level control.
3. **Fill LED** - Indicates when tank is filling.
4. **Pump Switch** - Activates dual speed circulation pump.
5. **Cool Switch** - Activates cooling temperature control.
6. **Cool LED** - Indicates temperature control is asking for cooling.
7. **Digital Display** - Displays water temperature and temperature control setting.
8. **Press-To-Set Button** - Switches digital display from reading the actual water temperature to reading the temperature control setting.
9. **Adjust Slot** - Adjusts (use small screwdriver supplied) temperature control settings when Press-To-Set button is pressed.

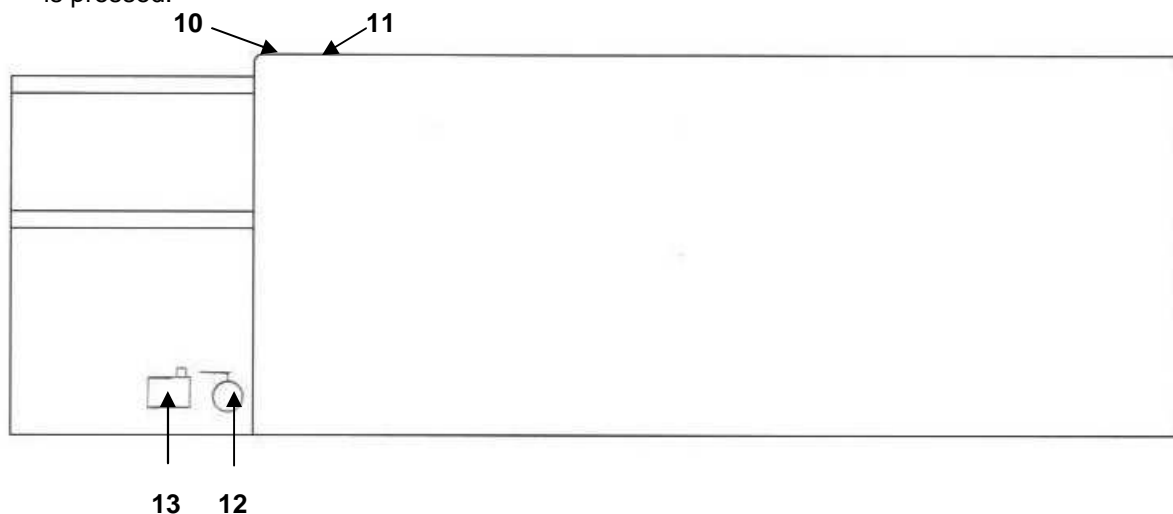

10. **High Speed Button** - Switches pump speed from low to high for 20 minutes.
11. **Air Control** - Controls air supply to hydrotherapy jets while pumps are running.
12. **Drain Valve** - Manually drains tank. Power switch should be switched off.
13. **Reset Button** - Resets refrigeration compressor high pressure cut-out following loss of cold water supply pressure.

**Installation:****1. Site Preparation**

The attached Spec Sheet indicates the utilities necessary for installation of the CRYOTherm. Consideration should also be made for the 4,685-pound load this unit will place on the floor and building structure.

**2. Level CRYOTherm**

Stainless steel hex head leveling feet are provided to give uniform support. Thread these into the frame before positioning the unit. Place a level on the top Corian surface and level the CRYOTherm.

**3. Drain Connections**

A 2-inch copper drain connection is provided. This should be positioned over an open, trapped floor drain with an adequate air gap. Refer to local plumbing codes for specific requirements.

**4. Cold Water Supply Connection**

A ¾" NPT connection is provided for connection to the cold water supply. An internal RPZ valve protects the supply and requires annual inspection. This supply fills the tank and cools the condensing unit.

When the supply pressure fails, the high-pressure cutout on the refrigeration system will trip making it necessary to press the reset button (see #13 under control functions) after the pressure is restored.

**5. Electrical Connection**

Power requirements for the Model 4902 CRYOTherm are 208-240/120 volt, single-phase, 60 hertz, 23 amps.

Only plug the CRYOTherm into a voltage-matched NEMA 14-30R 30 amp, 3-pole, 4-wire grounding outlet that is properly grounded, polarized, and ground-fault-circuit-interrupter (GFCI) protected.

Use professional electrical assistance and follow the National Electrical and local codes as a minimum.

## Operation:

### 1. Operator Considerations

Operators should have a working knowledge of clinical hydrotherapy and cryotherapy procedures. They should also have a complete understanding of the procedures described in this manual, the ability to operate the CRYOTherm, and assist patients.

### 2. Fill Tanks

If prescribed, add an antibacterial agent to the water at this time. Follow the antibacterial agent manufacturer's directions for use.

With Fill, Pump, and Cool switches "Off" on both tanks and drain valves "Closed" press Power Switch to "On". Now press Fill switches to "On" and allow water level to rise. An automatic level switch located behind the overflow grille will control the water level just below the overflow pipe height. These level controls have a ten-second delay with the under-level condition to avoid valve chatter.

### 3. Switch Pump "On" and Adjust Hydrotherapy Jets

When tanks are full, turn Air Control Valve to "Min" and press Pump switch to "On". Pumps will begin circulating the water. Pump must be left on at all times for proper operation of the temperature control. ***Without circulation ice will form on the tank walls.***

Tank has six adjustable hydrotherapy jets. When the dual-speed pump is switched "On" it circulates water through these jets in the low-speed mode. Press the white button to switch for high-speed hydrotherapy. After 20-minutes the pumps will resume low-speed operation automatically. Turn Air Control valve to "Max" for more vigorous aerated hydrotherapy.

### 4. Select Temperature

The temperature is adjustable from 40°F to 70°F. To set the temperature press the black Press-To-Set button and the digital display will switch from reading the water temperature to indicating the set point temperature. Using the small screwdriver supplied, adjust the set point to the desired temperature. When the button is released, the digital display will immediately return to reading the actual water temperature. Press the Cool switch to "On" and the CRYOTherm will seek the set point temperature and control at that point.

### 5. Patient Loading

CRYOTherm is ergonomically designed to accommodate all size patients and support them either standing in the tanks or sitting anywhere on the Corian surfaces.

**Maintenance:****1. Clean Regularly**

After patient use, clean the tanks thoroughly with a disinfectant and inspect for contaminants, particularly after treating open lesions.

**Note...**The entire CRYOTherm should be inspected, cleaned, and disinfected at least once a week.

**Caution...**When cleaning CRYOTherm wear protective gloves.

**2. Disinfection Procedure:**

- a) Turn Power switch "Off" and drain tanks before disinfection.
- b) Apply solution of diluted bleach (1/4 cup bleach to one gallon water) or EPA approved intermediate level disinfectant with cloth or sponge. Wet surface thoroughly and allow to remain wet for ten minutes.
- c) Rinse with clean water to reduce residue of disinfectant on the surfaces to a non-toxic level.
- d) Examine tank walls, exterior surfaces, and associated components for damage and visual cleanliness.

**3. Culturing**

Routine culturing of the tanks should be done. Guidelines for culturing can be obtained through clinical laboratory or hospital infection control department.

**4. Replace Filters:**

Filters should be changed following a positive culture or when circulation is impaired.

- a) Turn Power switch "Off" and drain tanks before replacing filters.
- b) Loosen locking ring (pull down locking tab and rotate ring counter clockwise)
- c) Remove old filter cartridge and clean inside of housing.
- d) Install new 100 sq. in. filter cartridge (Grimm part #060733) into filter housing.
- e) Clean o-ring. For ease of future disassembly, apply a very light film of silicone lubricant to o-ring and threads.
- f) Be sure o-ring is in place then and insert housing into filter head. Tighten locking ring (pull down locking tab and rotate ring clockwise)

**Service:****1. Trouble Shooting**

Problem	Possible Cause
Nothing happens when Power switch is switched "On".	Power cord unplugged. GFCI breaker tripped.
CRYOTherm won't cool.	Check refrigeration compressor high pressure cut out reset button.
Ice forms on tank walls.	Pump switch turned "Off".
Water circulation impaired.	Check pump suction fittings for debris. Check filter cartridges.

**2. One Year Warranty**

The CRYOTherm warranty applies to purchaser and any subsequent owner. Grimm Scientific Industries, Inc. (warrantor) warrants this product to be free from defects in material and workmanship under normal use for one year from the date the product is shipped. The warrantor may opt to replace the product, rather than repair it. This warranty does not cover damage resulting from accidents, improper use, or maintenance of this product, and in no case shall warrantor's liability exceed the original purchase price paid for the product. This warranty is expressly in lieu of all other express warranties.

**3. Repairs**

Should repairs be needed on your CRYOTherm, simply call the factory and ask to speak with a service technician. We will repair it free during the first year and quote a nominal parts and labor fee after the warranty expires.

**4. Inspection & Calibration**

Every year a qualified service technician should check the CRYOTherm. Temperature control calibration should be verified with a certified mercury thermometer, electrical system tested for grounding and leakage current, and RPZ valve inspected. Grimm can provide this service.

**5. Call 1-800-223-5395**

All Grimm products are sold and serviced directly by the factory. Call us toll-free Monday through Friday, 8:00 a.m. to 5: p.m. EDT about your thermal therapy needs. We are often available before and after hours so call us anytime.